



B350VT Series

URISAFE® 24-HR URINE COLLECTION CONTAINER
WITH DEVICE FOR VACUUM TUBES

- Sample volumes can be read in 100 ml increments



B350-3VT

B350-4VT



Contributing to the evolution of Disposable Laboratory Plasticware for over 47 years.

simport.com

B350VT Series

URISAFE® 24-HR URINE COLLECTION CONTAINER WITH DEVICE FOR VACUUM TUBES

Container made of high-density polyethylene, cap made of polypropylene



The innovative cap:

- Designed to use with vacuum tubes. It allows a clean transfer of the urine from the container to the vacuum tube.
- Snap valve pouring spout incorporated
- Leakproof: a cap liner ensures safe sample transport
- It includes a cannula with a needle inside protected by a rubber cap and a plastic suction nozzle. The cap features a label that seals the cannula in order to prevent contact with the collection needle.

The unique snap valve pour spout:


- Easily popped open or pushed shut
- Offers dripless pouring
- Controlled flow rate for better handling of poured volumes
- Reduced risks of aerosol contamination when pouring
- Eliminates splashing and exposure to hazardous body fluids

Large central handle:

- Can be gripped comfortably with three fingers
- Permits pouring off of samples with ease and reduces fatigue

Graduated vertically and horizontally:

- Graduations are easy to read
- Place container upright and read in 100 ml increments
- Can be gamma sterilized

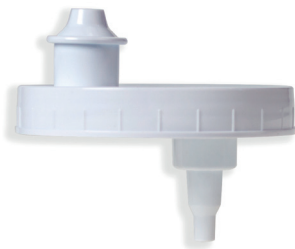

Cat #	Volume (Liters)	Dimensions (cm)	Qty/Cs
B350-3VT	3	11.5x24.5x16.0H	40
B350-4VT	4	11.5x 24.5x 20.1H	30